



Salsa® M



Salsa[®] M

Excellent indoor manoeuvrability with great outdoor performance

Excellent indoor manoeuvrability

- Quickie mid wheel drive technology delivers an ultra-compact turning radius for superior indoor manoeuvrability
- Narrow width of 60cm (23.6") enables easy access through doors and confined spaces
- Very low seat height (from 40cm (15.7"), and from only 43cm (17") with powered tilt) makes access under tables and into vehicles easy


Great outdoor performance

- Patented Quickie all-wheel independent suspension ensures a smooth and stable ride over rough or uneven terrain
- Quickie anti-pitch technology eliminates pitching on slopes, ensuring your safety at all times
- Effortlessly climbs up to a 7cm (2.8") kerb
- 50Ah gel batteries and powerful high-torque motors provides up to 26km (16 miles) of driving range at speeds of 6 or 10kph (4 or 6mph)



Fully adjustable seating system

The Salsa M seating system with width and depth adjustment ensures a perfect fit today and tomorrow.

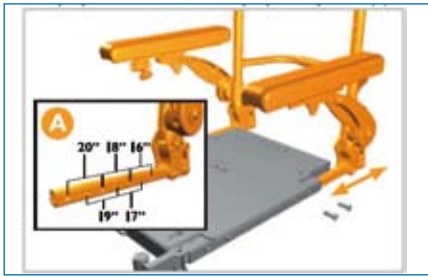

- Seat width and depth adjustment from 41 to 51cm (16 - 20")
- Backrest height adjustable in three positions from 45.5 to 50.5cm (18 - 20")
- Five adjustable seat heights and four seat angles
- Easy adaptation for alternative seats, backs and cushions (e.g. standard back sling, Jay, third party seating)



NEW: Low push force, emergency freewheel allows the Salsa to be moved very easily when it is in freewheel



Flip-up armrests and foot rests for easy accessibility



Lift-off backrest for easy transportation and storage

Optional Jay Comfort Seating System

A variety of Comfort Seating options are available, these include:

- Comfort seat pan
- Standard and optimum comfort cushion
- Shallow and deep contour backrests
- Comfort calf pads
- Comfort and flexible headrests



Optimum comfort cushion



Deep contour backrest



Comfort headrests



Advanced electronics

VR-2 electronics as standard

- 70 Ah controller as standard. 90 Ah heavy duty version optional
- Dual control (option)

R-net expandable controls (optional)

- Provides simple programming options
- Shortcut to your favourite functions at the touch of a button thanks to Quickie's patented assignable button technology
- Control up to 6 powered functions through your joystick
- Drive profiles can be tailored to your personal requirements
- Control your PC with your joystick via Bluetooth® (option)
- Dual control (option)

Full range of powered options available

Available powered options:

- C.O.G tilt-in-space (0° to 30°)
- Lift and tilt (30cm (11.8") lift, 0-30° tilt)
- Recline
- ELR/ALR powered leg rests
- Easy adaptation of powered options if required (upgrade kit)

Optional retro-fit kits:

A full range of upgrade kits can be added to your chair at a later date if your needs change

Salsa[®] M

Freedom indoor and out



Technical Data: Seating

Seat angle:	Adj. 0°, 3°, 6°, 9°	
Seat depth:	41, 43, 46, 48, 51cm (Adjustable from 41 to 51cm)	16, 17, 18, 19, 20"
Seat width:	41, 43, 46, 48, 51cm (Adjustable from 41 to 51cm)	16, 17, 18, 19, 20"
Seat height:	40.5 - 45.5cm (no tilt), 45.5 - 47.5cm (powered tilt), 48cm (lift & tilt)	16 - 18" (no tilt), 18 - 18.5 (powered tilt). 19" (lit & tilt)
Backrest angle:	Adj. -4°, 0°, 4°, 8°, 12°, 16°	
Backrest height:	Adj. 46, 48 and 51cm	18, 19 and 20"
Backrest width:	Adjustable from 41 to 51cm	16 - 20"
Lower leg length:	33 - 51cm	13 - 20"
Legrest angle:	70°, 90°, ALR, ELR	
Armrest height:	Adj. from 22 to 31cm	8.7 to 12"
Armrest depth:	Adj. from 31 to 41cm	12 to 16"

Specification

Overall width:	60cm (std wheels), 63cm (low push force wheels)	23.5" (std wheels), 24.5" (low push force wheels)
Overall length: (including legrests)	110cm	43"
Chair weight*:	113 - 165kg	249 - 363lbs
Turning circle:	120cm	47"
Battery capacity	2 x 50 Ah	
Range:	26 km (ISO 7176-4)	16 miles
Speed:	6kph and 10kph	3.7mph and 6mph
Kerb climbing ability:	7.5cm	3"
Maximum user weight:	140kg	22 stone
Max safe slope:	8°	14%
Ground clearance:	5cm	2"
Transit tested to ISO 7176 Part 19 and passed:	Yes	
Warranty:	2 years on frame and electronics, 3 years on motors, 1 year on batteries (conditions apply) and chargers	

Find out more today!

Visit www.quickiepower.com for more information, to find your local dealer or to request a demonstration.



Get SOCIAL with QUICKIE

Keep up-to-date with the latest news on the Quickie Salsa M and other Quickie products through our Facebook page:

www.facebook.com/QuickieWheelchairsUK



Sunrise Medical Limited
Thorns Road
Brierley Hill
West Midlands
DYS 2LD

Tel: +44 (0) 845 605 66 88
Fax: +44 (0) 845 605 66 89
www.SunriseMedical.com

* Please note: Total weight of the powerchair depends on the configuration
113kg (base chair) to 165kg (base chair plus solid tyres, power lift, tilt, recline, headrest, comfort seating, R-net controls and lights)