




WHERE YOU WANT TO BE.

Sapphire[®] 2



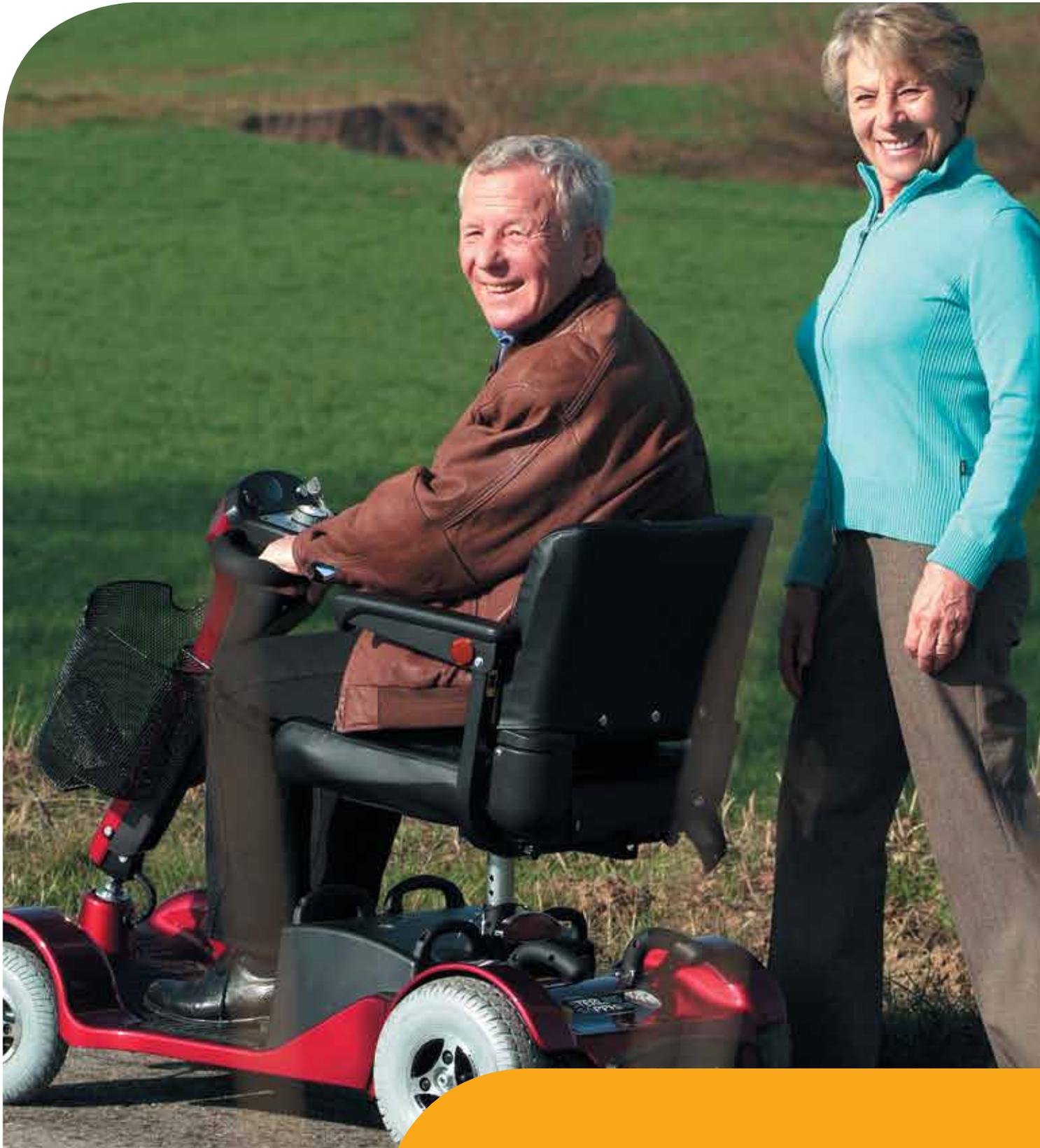

Sapphire[®] 2

Take it apart,
take it away,
take it further!





Detachable

To get your Sapphire 2 ready for transporting, simply lift off the seat, lift out the batteries, and lift the rear handle to separate the front and rear of the scooter. It really is that simple! By the way, it goes back together just as smoothly!



The Sterling Lock system is cable and connector free so there are no fiddly bits to deal with.

Transportable

Now the Sapphire 2 can be loaded into your vehicle. Luckily, as one of the lightest scooters in its class, this is no problem with the Sterling Lift dual handles fitted to the batteries and rear drive unit. These make lifting and handling the parts safer and more comfortable.



With its large batteries carrying you in style and comfort, the Sapphire 2 is the scooter you'll want to go further on!


The tiller can be safely secured in position for transporting.

Comfortable

Now you're ready! Stretch out your legs and take advantage of the Sterling Comfort package: three wheel legroom with four wheel stability, infinitely adjustable tiller, wrap-around handle bars and finger-thumb directional control all fitted as standard.



Sapphire[®] 2

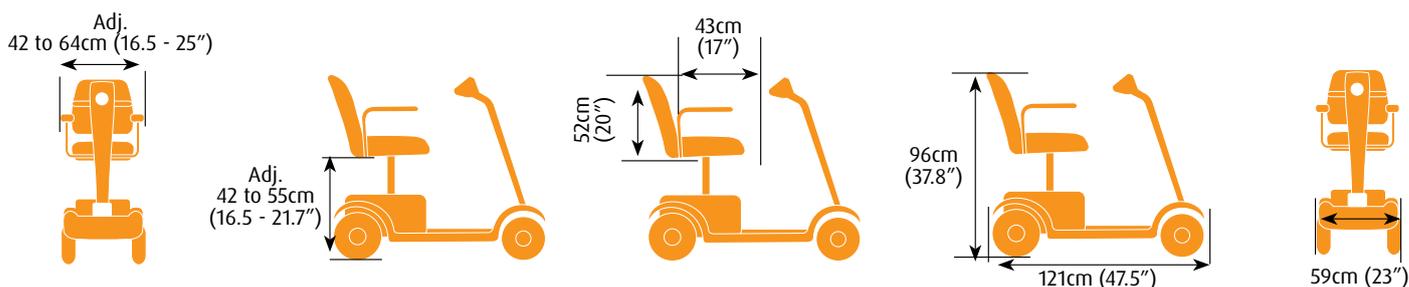

Technical data sheet

Version	4 Wheels	Drive unit weight	19.5kg (43lbs)
Speed	6.4kph (4mph)	Battery weight	unit 12kg (26.4lbs)
Batteries	35 Amp	Total weight when folded/dismantled	22.5kg (49.5lbs)
Charger off board	5 Amp	Turning radius	140cm (55")
Maximum range	30km (19 miles)	Maximum safe slope	14° - 113kg (17 stone) user weight / 10° - 150kg (23.8 stone) user weight
Colour	Blue & Red	Ground clearance	9.5cm (3.7")
Maximum user weight	150kg (23 stone)	Tyres	Pneumatic
Total weight	79kg (174lbs)	Front wheel size	22.5cm (9")
Heaviest part (front section)	22.5kg (49.6lbs)	Rear wheel size	26.5cm (10.5")
Seat weight	12.5kg (27.5lbs)	Stylish lights	Yes

Find out more today!

Visit www.sterlingscooters.co.uk for more information, to find your local dealer or to request a demonstration.

For further information on the full specification, options and accessories, please refer to the order form. All information is subject to change without notification. Please consult Sunrise Medical with any queries you may have.

If you are visually impaired, this document can be viewed in PDF format at www.SunriseMedical.co.uk or alternatively is available on request in large text.



Sunrise Medical Limited
Thorns Road
Brierley Hill
West Midlands
DY5 2LD
Tel: +44 (0) 845 605 66 88
Fax: +44 (0) 845 605 66 89
www.SunriseMedical.co.uk